

# Deflection of Electrons

Every statement in physics has to state relations between observable quantities.

E. Mach (1838-1916)

## OBJECTIVES

To determine the effect of electric and magnetic fields on a beam of electrons. To measure the charge to mass ratio of an electron.

## THEORY

In this exercise we ask how a stream of negatively charged particles responds to electric and magnetic fields, and thereby determine the speed of the particles and their charge to mass ratio. Consider a beam of particles with charge  $e$  and mass  $m$ , moving at speed  $u$  through a vacuum. If there is a uniform electric field  $\vec{E}$  and uniform magnetic field  $\vec{B}$ , The force on each particle is given by

$$\vec{F} = e\vec{E} + e\vec{u} \times \vec{B} \quad (1)$$

We might choose to produce the  $\vec{E}$  field with a voltage applied to two parallel plates, and to produce the  $\vec{B}$  with a coil outside the vacuum chamber. In either case, the fields are calculable from the geometry, shown in Fig. 1, and measurable voltages and currents. For simplicity, we consider only one field at a time, and calculate the deflection of the beam when we turn on each field. In this analysis, we are neglecting the effect of gravity since the acceleration due to the electric and magnetic field used in this lab will be much times greater than acceleration due to gravity.

The electric field causes a constant transverse acceleration during the time the particles are between the plates. The force acting is simply  $eE$ , and the time is  $\Delta t = L_p/u$ , so the particles acquire a velocity  $u_1$ , perpendicular to  $u$ , given by

$$u_1 = \frac{F}{m} \Delta t = \frac{e}{m} \frac{EL_p}{u} \quad (2)$$

After leaving the plates the particles travel in a straight line to the end wall of the tube. The path makes an angle  $\theta$  with the original direction, resulting in a displacement  $D_E$  given by

$$D_E = L_E \tan \theta = L_E \frac{u_1}{u} \quad (3)$$

The electric field can be expressed in terms of the voltage  $V_d$  applied to the plates by  $E = V_d/a$ .

Using this result in Eq. 2 and 3, we get

$$D_E = \frac{e}{m} \frac{L_E L_P}{au^2} V_d \quad (4)$$

Since the tube wall glows where the beam hits, we can measure  $D_E$  with a ruler. (There is also a small displacement that occurs within the plates, given by  $eE(\Delta t)^2/2m$ . This is only about 5% of  $D_E$  for the geometry we will use, and is neglected.)

From Eq. 1, the magnetic force is constant and perpendicular to  $\vec{u}$ . The particles will move in a circle whose radius may be found by equating the magnetic force to the centripetal force needed to move in a circle of radius  $R$

$$euB = \frac{mu^2}{R} \quad (5)$$

or

$$R = \frac{mu}{eB} \quad (6)$$

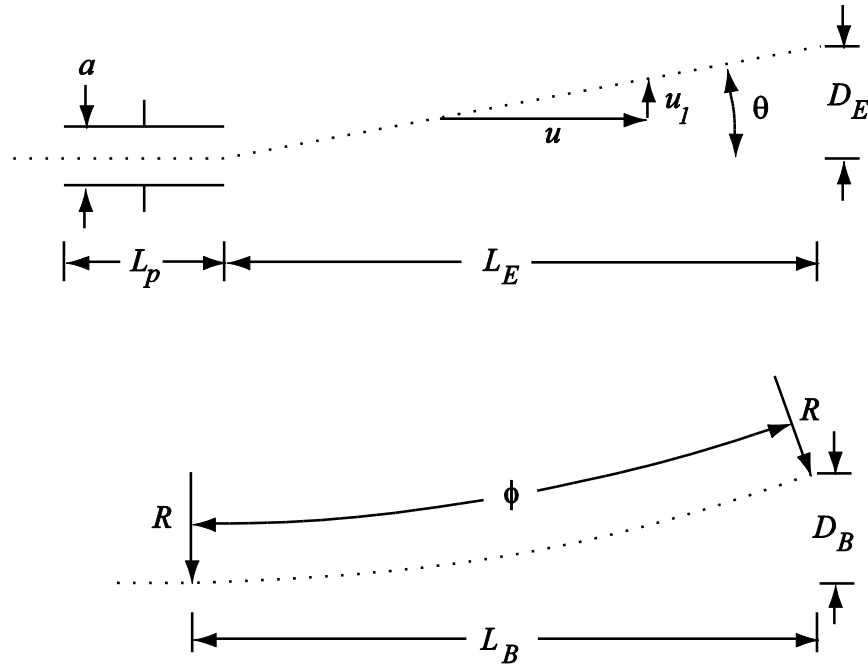


Fig 1 Arrangements for electric or magnetic deflection of charged particles. The magnetic field is uniform over the entire distance  $L_B$ , while the electric field exists only between the plates.

Referring to Fig. 1 again, the total deflection  $D_B$  is given by

$$D_B = R - R \cos \phi \quad (7)$$

For small angles, we can use the approximations  $L_B \approx R\phi$  and  $\cos \phi \approx 1 - (\phi^2/2)$  to combine Eq. 6 and Eq. 7 into

$$D_B = \frac{e}{m} \frac{L_B^2}{2u} B = \frac{e}{m} \frac{L_B^2}{2u} GI \quad (8)$$

In the last part we expressed the field  $B$  by  $B = GI$ , where  $I$  is the current in the coils and  $G$  is a geometric factor that will be given later.

Eqs. 4 and 8 are our main results. **They tell us that graphs of  $D_E$  vs  $V_d$  and  $D_B$  vs  $I$  will be straight lines with specified slopes. Calling the slopes  $s_E$  and  $s_B$  respectively, we can solve for the two unknowns  $u$  and  $e/m$ , to arrive at**

$$u = \frac{2L_p L_E s_B}{a G L_B^2 s_E} \quad \text{and} \quad \frac{e}{m} = \frac{4L_p L_E s_B^2}{a G^2 L_B^4 s_E} \quad (9)$$

These two equations will be used for analyzing our data.

## EXPERIMENTAL PROCEDURE

Figure 2 is a sketch of the vacuum tube that we will use, showing the electrical connections needed. Electrons are boiled out of the cathode with a heater. Various electrodes, labeled as grids and anodes, form the stream of electrons into a small focused beam. **The voltage labeled C- is adjusted to provide the proper confinement for the electron beam, while the B+ and C- voltages together determine the energy and speed of the electrons.** After leaving the accelerating and confining electrodes, the beam passes between two sets of deflecting plates which are arranged to deflect the beam in orthogonal directions. As shown in Fig. 2, we will use only the second set, the set nearer to the screen. The other set is connected to the COM (i.e. common) terminal so it will not charge up and produce undesirable deflections. After the deflection plates, the beam travels a distance  $L_E$  to the end of the tube, where it causes a special material to **glow green** at the point of impact. The coils which produce the magnetic field will be discussed later. Table 1 near the end of the manual contains a summary of the constants you will need for analyzing your results.

Before proceeding with the experiment, there is a hazard you must be aware of. The power supply used here produces dangerously high voltage. Because of this it has two control switches. The AC POWER switch turns on the power supply and the low voltage outputs. It may be left on throughout the experiment. **The DC ON/STANDBY switch controls the high voltage.** **This switch must be at STANDBY when you make any changes in the wiring. Failure to follow this procedure may result in a dangerous electrical shock.** It is good practice to set the switch to STANDBY any time you are not actually making measurements.

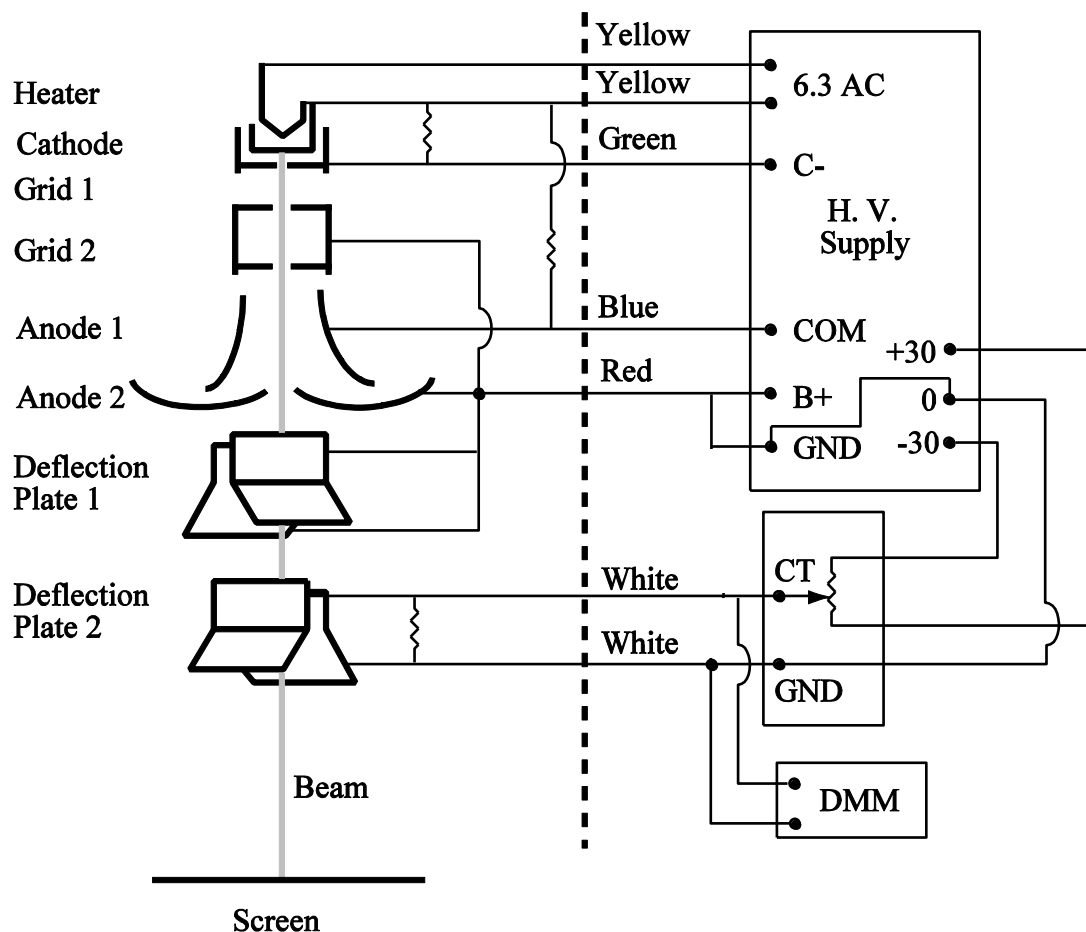


Fig. 2. Electrode configuration and wiring diagram for the cathode ray tube (CRT). **Most of the wirings between the CRT and the HV Supply should already be in place.** The only wiring you should have to do is to connect the white wires to the small box with a knob marked with CT and GND, which contains the circuit driving the deflection plates, and to connect the DMM in parallel with the white wires to measure the deflection voltage. **Wires to the DMM are connect to the terminals marked VΩ and COM, and to measure the deflection voltage, set the large knob in the center of the DMM to point to the V with the solid and dashed lines over it.**

Check that the tube and deflection circuits are connected as shown in Fig. 2. **Follow the wiring instruction given in the caption to figure 2.** Turn on the AC POWER switch and set the deflection voltage,  $V_d$ , to close zero. **You adjust the deflection voltage,  $V_d$ , by turning the knob on the small box that is marked with CT and GND, after you had connected the white wires from the CRT to the small box.** Allow the heater to warm up for about a minute, and then turn on the DC voltage with the DC ON switch. The VOLT METER switch on the H.V. Supply connects the meter on the power supply to either the **B+ VOLTS (red scale)** or the **C- VOLTS (black scale)**. Note that the B+ and C- potentials are relative to the COM (i.e. common) terminal. This means that, for example, anode 1 will be at a potential B+ more negative than anode 2 (which is at ground potential). First, set the B+ potential to about 250 V, and then adjust the C- potential to get the sharpest possible (smallest) spot on the screen of the tube. **Next, try several settings of the B+ and C- to find the largest value of B+ for which you can get the smallest spot, and then try to find the lowest value of B+ for which the spot is small and bright enough to easily see. You will have to adjust the C- voltage to get the smallest spot at each B+ setting.** Note the largest and smallest values of B+ that your tube allows. These correspond to the range of speeds  $u$  that are available, and will be used in the deflection measurements. **A larger value of B+ corresponds to an electron beam with faster electrons.**

There are two extraneous effects that you may notice. The first is that the position of the beam spot depends on the orientation of the tube. This occurs because the earth's magnetic field can deflect the beam. Pick an orientation which places the spot conveniently and do not move the tube during this part of the experiment. **The second effect is more common and it is the appearance of a large black area on the screen when operating at lower B+ voltages or when the electron beam is relatively large. At higher B+ voltages, electrons will be ejected from the screen as fast as they arrive, maintaining charge neutrality on the screen. At lower B+ voltages or when the electron beam is relatively large, the screen may become charged, thereby defocusing and repelling the beam. The problem can be cured by temporarily increasing the B+ voltage to the maximum value then sweeping the deflection voltage over the available range several times. This will cause the electron beam to move vertically on the screen. You can adjust the deflection voltage by turning the knob on the small box with the knob connected to the apparatus via the white wires. This will remove the charges on the screen. Be sure to lower the B+ voltage to your desired value after the dark area has been removed.**

### *Electric deflection*

With the **largest value of B+ voltage for which you can get the smallest spot**, increase the deflection voltage  $V_d$  by turning the knob on the small box. The spot on the tube face should

move vertically in response to the electric field. Measure the deflection  $D_E$  for **at least five values of  $V_d$** , and plot  $D_E$  as a function of  $V_d$ . The values of  $D_E$  corresponds to the vertical distance traveled by the green spot on the screen relative to a fixed position on the screen of your own choosing. **You can measure  $D_E$  with a ruler or the equally spaced lines on the screen. The values of  $V_d$  are measured using the DMM that you have connected to the small box.** Reset the B+ voltage to the lowest value that you found earlier (Be sure to adjust the C- voltage to get the smallest green spot) and repeat the measurement, plotting the new data on the same graph. The slopes of these two lines are the values of  $s_E$  for two different values of  $u$ . Set the power supply to STANDBY and proceed to the magnetic deflection. **Which setting of B+ gives you the larger value of  $s_E$ ? Explain the physical reason for your answer.**

### *Magnetic deflection*

To prevent deflection plate 2 from accumulating a negative charge, **connect both white wires from the tube that were connected to the small box with the knob to the GND terminal on the high voltage power supply.** This procedure assures that there will be no stray electric fields to deflect the beam.

Position the tube on the plywood platform between the coils so that **the region between the end of the electron gun and the screen is centered between the coils.** The magnetic field will be approximately uniform in this region. **Connect the coils, the low voltage power supply (This is the same power supply with a single knob that you used for the two pervious labs), and the DMM in series so that current flows sequentially through each. The DMM will measure current flowing through the terminals labeled mA and COM when the large knob at the center of the DMM points at the position labeled mA.** The field is proportional to the current read on the DMM, with constant of proportionality  $G$  listed in Table 1.

|       |                                |
|-------|--------------------------------|
| $L_p$ | $1.5 \times 10^{-2}$ meter     |
| $L_B$ | $1.6 \times 10^{-1}$ meter     |
| $L_E$ | $1.2 \times 10^{-1}$ meter     |
| $a$   | $4.0 \times 10^{-3}$ meter     |
| $G$   | $4.1 \times 10^{-3}$ Tesla/Amp |

Table 1 Tube and Coil Constants

Turn on the high voltage DC, with no current in the coils, and set the B+ to the higher value you used before (Be sure to adjust the C- voltage to get the smallest green spot). As you turn up the current in the coils by turning the knob on the low voltage power supply, the spot

should deflect. If the spot moves only slightly or becomes defocused, reverse the connections on one of the coils and try again. Measure the deflection  $D_B$  for **at least five values of current  $I$** , and plot  $D_B$  as a function of  $I$ . Repeat the procedure for the lower B+ voltage you used before (Be sure to adjust the C- voltage to get the smallest green spot), plotting the data on the same graph. When you have both sets of data, turn off the high voltage with the DC ON/STANDBY switch.

### *Analysis*

You now have values of  $s_E$  and  $s_B$  for two different values of  $u$ . Use Eq. 9 and the constants in Table 1 to obtain **two estimates of  $e/m$  and two values of  $u$**  from your data. If you use SI units throughout, your answers will automatically come out in coulombs/kilogram and meters/sec, respectively. For the conditions we are using,  $u$  will be about  $10^7$  m/s and of course will increase with increasing B+ voltage. The accepted value of  $e/m$  is  $1.76 \times 10^{11}$  coulombs/kg, and should be the same for any value of  $u$ . Your number may not be very precise but it will be clearly different than  $e/M = 9.42 \times 10^7$  coulombs/kg found for hydrogen in chemical experiments.